
COAGULATION



MOZAÏK™

Microperforated adhesive tape post-dialysis compression device for guaranteed hemostasis



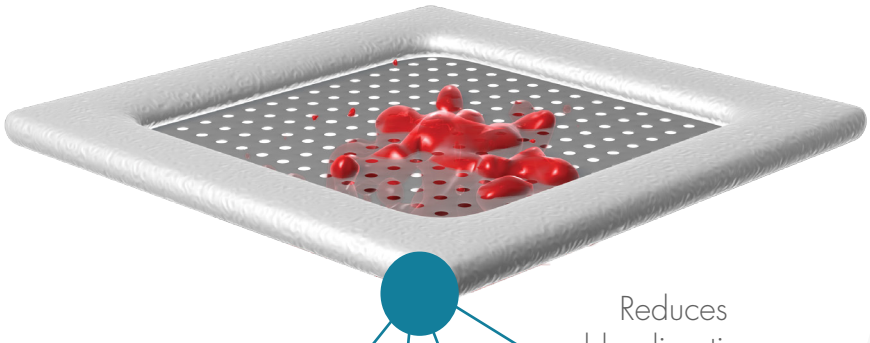

- Reduces bleeding time
- Enables real-time monitoring of the coagulation process
- Thick foam frame for easy installation and prevention of blood leakage
- Free of chemical agents and hemostatic substances



-54% **Reduces compression time by 54%**
*vs standard protocols**

MADE IN FRANCE

Mozaïk™ in reality

Discover now the Mozaïk™ device in video

nephrokit innovation

*According to EVAMOZ study, ECHO, 2024

NEPHROKIT, chez Gergonne Industrie
ZI nord - Rue de Tamas - CS 70204 01117 Oyonnax Cedex - France
www.nephrokit.fr


Transform your hemostasis process with our sterile, single-use Mozaïk™ device. Equipped with a patented microperforated adhesive film, Mozaïk™ acts blood rheology to fast-track platelet plug formation and cut compression time down to an average of just 5.8 minutes*, delivering unrivaled efficiency backed by scientific proof.
With Mozaïk™, safety and transparency go hand in hand. Our protective film shields against blood exposure accidents while allowing you to keep a close eye on the coagulation process. The thick foam edges make Mozaïk™ setup a breeze, guarantee a perfect seal, and streamline needle removal – ideal for self-dialysis or home dialysis.
Designed for compatibility with all Arterio-Venous Fistula (AVF) types, Mozaïk™ is committed to improving patients' quality of life. Our biocompatible adhesive is free from chemical agents, ensuring complete respect for patients' skin and health.
Opt for Mozaïk™, the innovative and efficient solution for rapid, safe, and comfortable hemostasis in dialysis.

Find out more


**According to EVAMOZ study, ECHO, 2024*

Contact us for further information


HOW TO USE IT

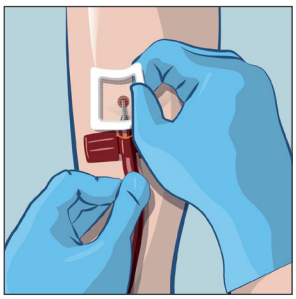


Apply Mozaïk™ to the puncture site after partially removing the needle.

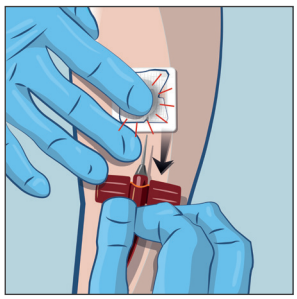

Place a sterile absorbing gauze on top of the Mozaïk™ and remove the needle entirely.

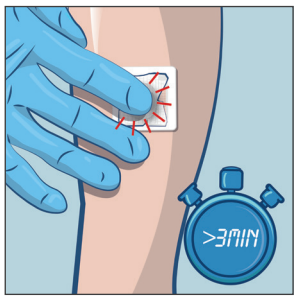

Compression for 3 minutes and control of the bleeding point.

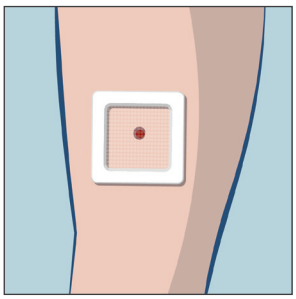

Assess hemostasis every minute. When coagulation is achieved, either use Mozaïk™ as the end-of-dialysis dressing or add another dressing as needed.

Mozaïk™ can be removed after 2 hours and within 24 hours at most.

DISTRIBUTED BY

HEMATRONIC

HEMATRONIC SRL
Str. Câmpului 253/2, Cluj-Napoca, Romania
Tel. +40 744 987 316
office@hematronic.ro
www.hematronic.ro

07/2024

nephrokit
innovation

NEPHROKIT, chez Gergonne Industrie
ZI nord - Rue de Tamas - CS 70204 01117 Oyonnax Cedex - France
www.nephrokit.fr