
FIXATION




PAPYRUS™

Adhesive medical device to reduce the risk of AVF needle dislodgement.



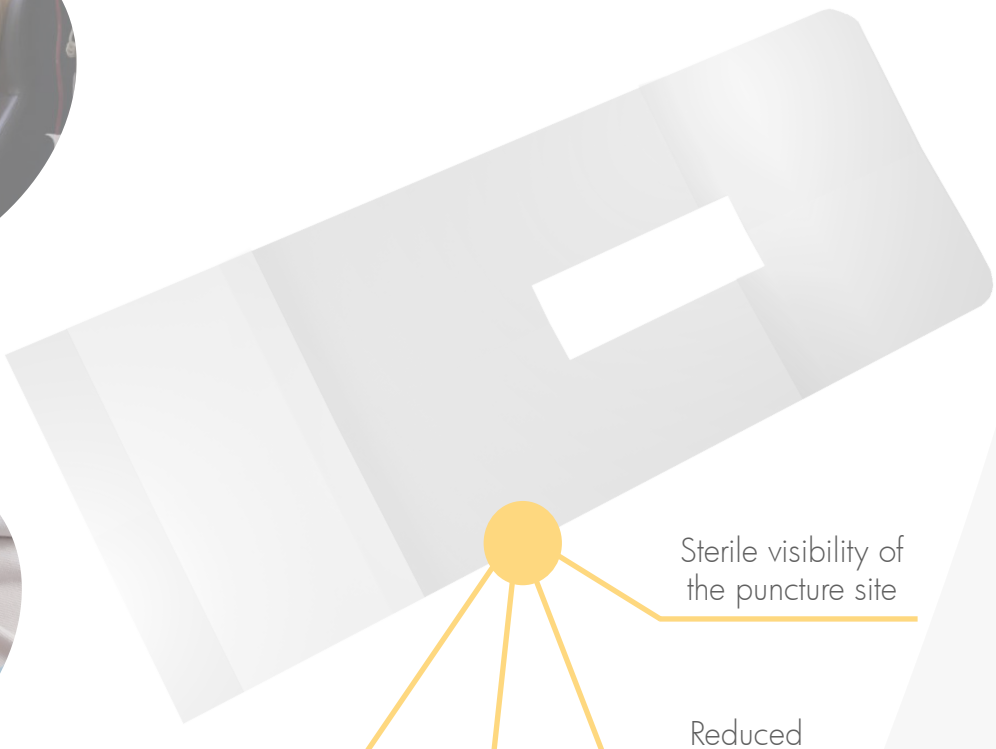

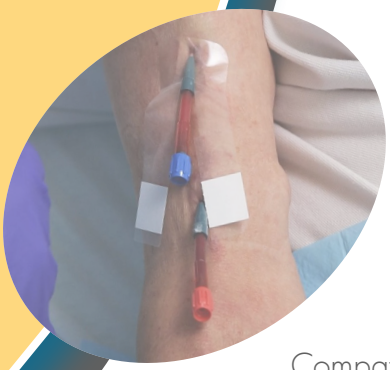

- Sterile visibility of the puncture site
- Reduced adhesive surface
- Compatible with all types of AVF and catheters
- Strong catheter fixation - prevents from dislodgement

Use with needle

Use with catheter

Papyrus™ in reality

Discover now the Papyrus™ device in video

NEPHROKIT, chez Gergonne Industrie
ZI nord - Rue de Tamas - CS 70204 01117 Oyonnax Cedex - France
www.nephrokit.fr

nephrokit
innovation


Developed by a Nephrologist, Papyrus™ is a sterile single-use device designed to reduce the risk of catheter dislodgement from arteriovenous fistulas. Universal, Papyrus™ is suitable for both needles and soft catheters, making it possible to standardize the installation protocols.

Easy to use, Papyrus™ makes the actions of nursing staff safer, faster and more efficient, although their habits are still respected, since fixation by tie remains possible.


Improving patients' quality of life is our priority, that is why Papyrus™ is made of a transparent film that ensures sterile occlusion of the puncture site without interfering with its monitoring. To limit the risk of irritation, the Papyrus™ adhesive is biocompatible and we reduced the surface in contact with the skin to a minimum.

HOW TO USE IT


• Use with fistula needle

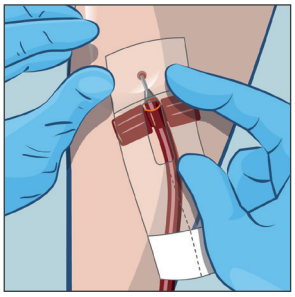


Fix the Papyrus™ on the puncture site using the adhesive featuring a notch.

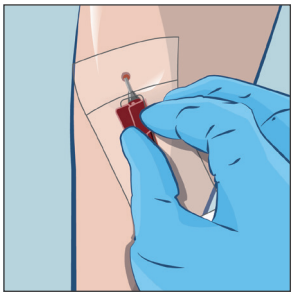

Fold the flexible wings and run them through the slot.

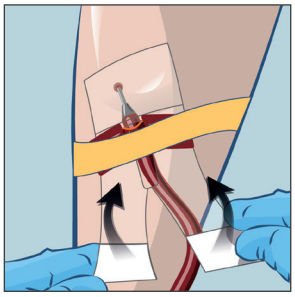


Place an adhesive strip over the wings or the base of the flexible catheter. Separate the 2 arms along the cutting line.

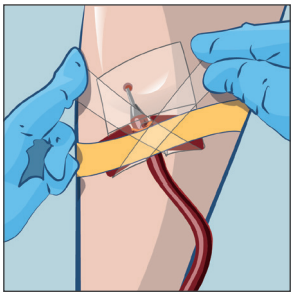


Fix the 2 arms on the skin by crossing over the wings.


• Use with fistula catheter

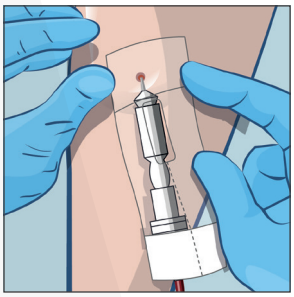


Fix the Papyrus™ on the puncture site using the adhesive featuring a notch.

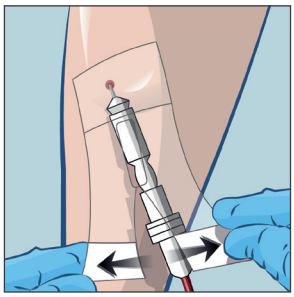

Separate the 2 arms along the cutting line.

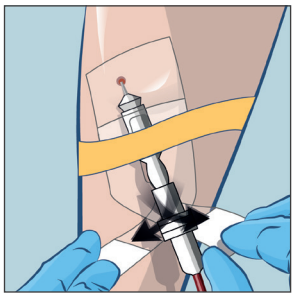

Put an adhesive strip over the fistula catheter.

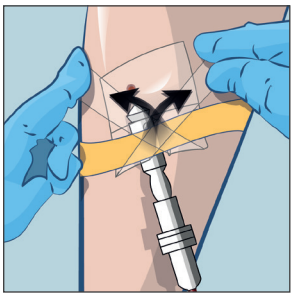

Cross the 2 arms over the base of the fistula catheter and fix them on the skin.

DISTRIBUTED BY

HEMATRONIC

HEMATRONIC SRL
Str. Câmpului 253/2, Cluj-Napoca, Romania
Tel. +40 744 987 316
office@hematronic.ro
www.hematronic.ro

11/2022

nephrokit
innovation

NEPHROKIT, chez Gergonne Industrie
ZI nord - Rue de Tamas - CS 70204 01117 Oyonnax Cedex - France
www.nephrokit.fr