PROTECTION








SolüKAT™

All-in-one medical adhesive device for Central Venous Catheter (CVC) protection.

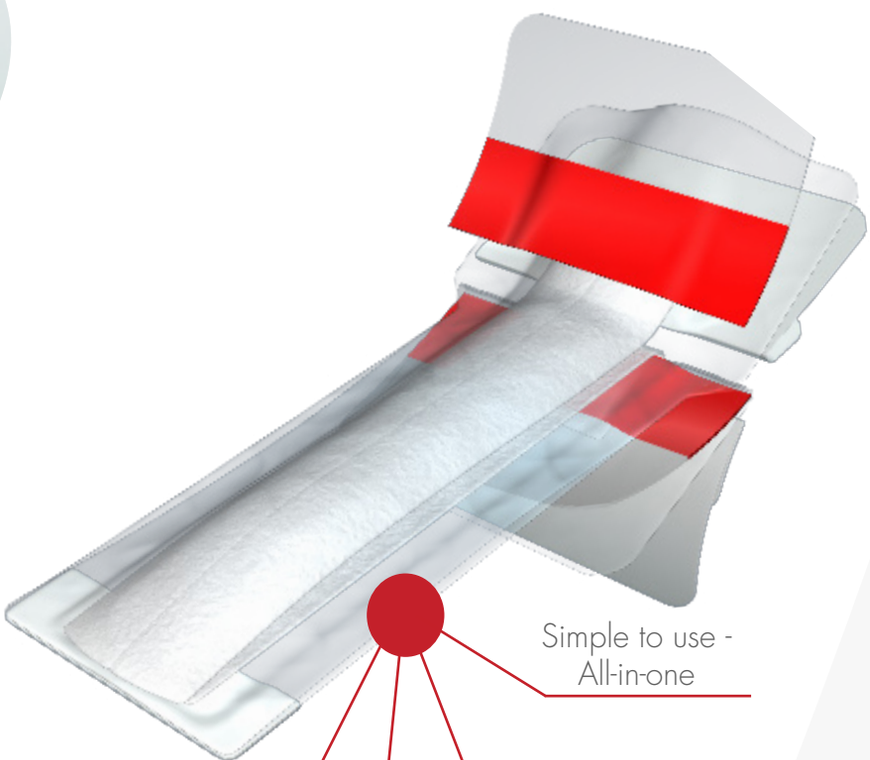

- Simple to use - All-in-one
- Very comfortable - soft and discreet
- Reduces the risk of exit-site infection
- Waterproof, allows showering

Solukat™ in reality

Discover now the Solukat™ device in video

MADE IN FRANCE

nephrokit innovation

NEPHROKIT, chez Gergonne Industrie
ZI nord - Rue de Tamas - CS 70204 01117 Oyonnax Cedex - France
www.nephrokit.fr

Solukat™ is an easy-to-use, single-use, sterile device designed to fully protect a central venous catheter in compliance with the strictest hygiene protocols.


Solukat™ is made of polyurethane, a material with a soft touch and unique characteristics since it is simultaneously a bacterial barrier, waterproof (according to the test of annex J of NF EN 13726:2023 standard) and highly permeable to gases (breathable).

Applied at the end of the dialysis session, Solukat™ can remain on the skin for up to 3 days, providing adequate protection from the exit sites to the connectors until the next session.

Improving patients' quality of life is our priority, that is why Solukat™ is a one-piece device perfectly sealed that allows the patient to enjoy showering again. To limit the risk of irritation, Solukat™ adhesive is biocompatible, and we also reduced the surface in contact with the skin to a minimum.


Patients with the most sensitive skin should use Solukat Si™, provided their skin does not perspire profusely, naturally, or as a result of sporting activities.

HOW TO USE IT

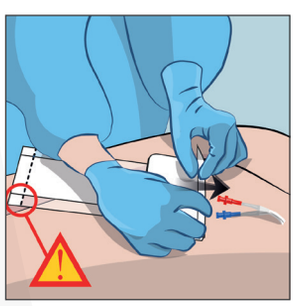


Slide the catheter ends into the bag and up to the insertion points.

⚠ Place the adhesive strip on the skin side.

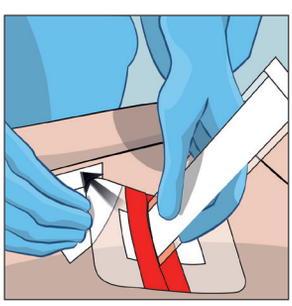


Remove the protective paper and stick the upper part.

Remove the protective paper from the lower part.

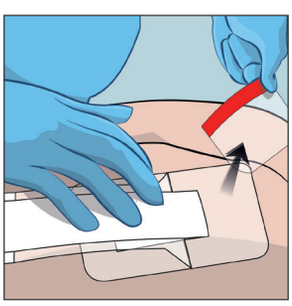

Remove the red marking strips.

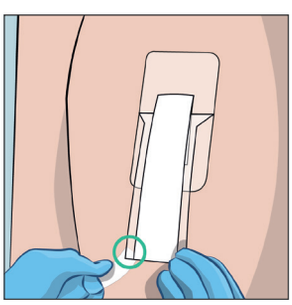

Remove the protective paper from the adhesive strip. Solukat™ is now ready for use.

DISTRIBUTED BY

HEMATRONIC

HEMATRONIC SRL

Str. Câmpului 253/2, Cluj-Napoca, Romania

Tel. +40 744 987 316

office@hematronic.ro

www.hematronic.ro

06/2026

NEPHROKIT, chez Gergonne Industrie

ZI nord - Rue de Tamas - CS 70204 01117 Oyonnax Cedex - France

www.nephrokit.fr

nephrokit innovation